

U+File

Safe Classic Endo




NIC[®]

For Better Endo


CAN) **ENT**
GROUP INC.

U+File

S a f e C l a s s i c E n d o

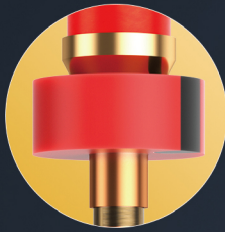


Radiopaque stopper




Precise depth marking



Continuous rotation

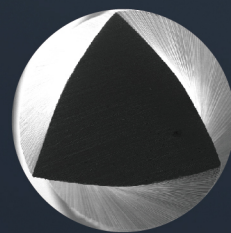

Convex triangle

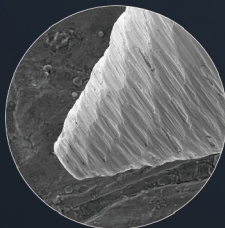


Safe guiding tip

Features

- MaxTech 3.0 NiTi alloy
- Triangular cross section design
- Variable taper design to have accurate shaping
- Perfect balance between flexibility and cutting ability
- Extended fatigue life
- Safe classic sequence


								
U+File	SX	S1	S2	F1	F2	F3	F4	F5
Size	19/.04	18/.02	20/.04	20/.07	25/.08	30/.09	40/.06	50/.05
Speed(RPM)	300	300	300	300	300	300	300	300
Torque(N.CM)	3.0	1.5	1.5	1.5	1.5	1.5	1.5	3.0
GP				F1	F2	F3	F4	F5

Sequence

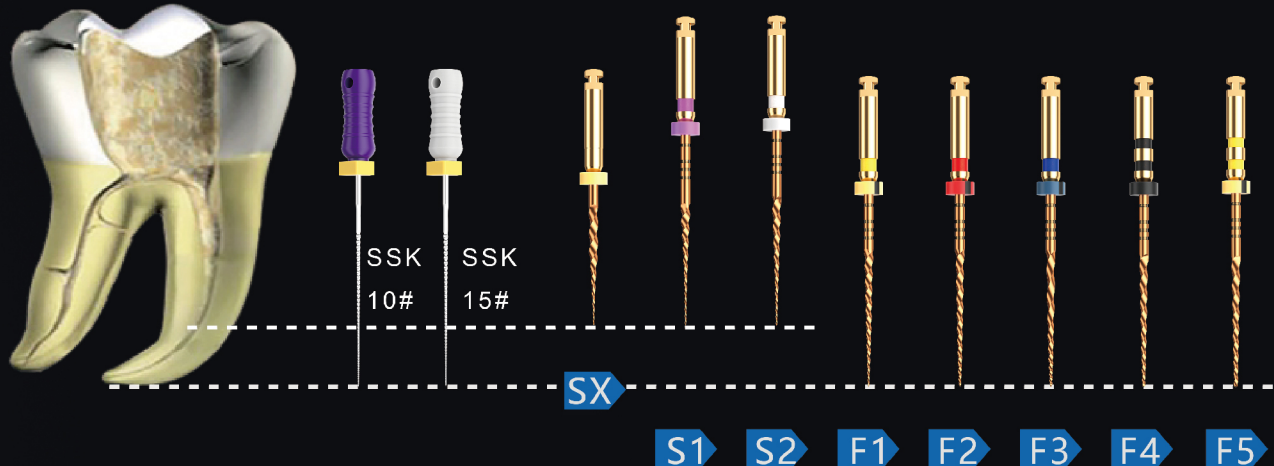



- 1 Scout the canal with SS K #10
- 2 Explore the canal with SS K #15
- 3 Use SX to expand the canal orifice
- 4 Scout coronal 2/3 & shaping this part with S1, S2
- 5 Scout apical 1/3 & follow S1, S2 until working length
- 6 Finish the apical 1/3 with F1, F2
- 7 F4, F5 followed when needed

References

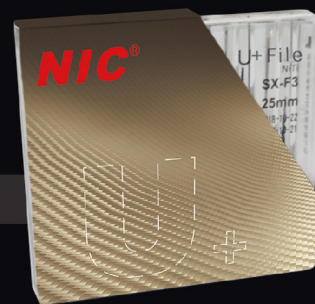

Size	19mm	21mm	25mm	31mm
SX	AFU 019 0SX			
S1		AFU 021 0S1	AFU 025 0S1	AFU 031 0S1
S2		AFU 021 0S2	AFU 025 0S2	AFU 031 0S2
F1		AFU 021 0F1	AFU 025 0F1	AFU 031 0F1
F2		AFU 021 0F2	AFU 025 0F2	AFU 031 0F2
F3		AFU 021 0F3	AFU 025 0F3	AFU 031 0F3
F4		AFU 021 0F4	AFU 025 0F4	AFU 031 0F4
F5		AFU 021 0F5	AFU 025 0F5	AFU 031 0F5
SX-F3		AFU 021 A01	AFU 025 A01	AFU 031 A01

Clinic Cases

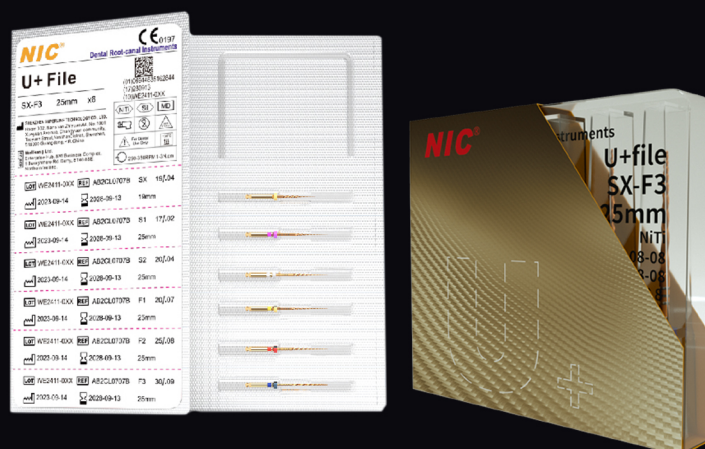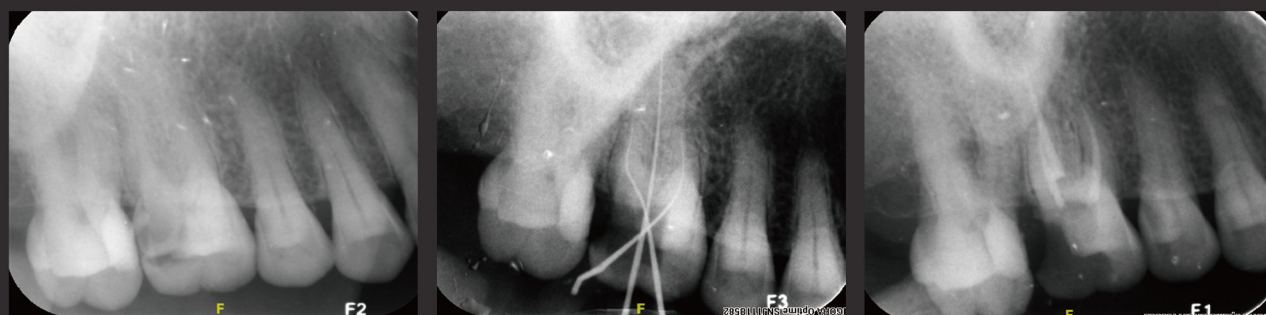

Available in sterile pack.



☎ 437-347-7791

✉ info@candentgroup.ca

🌐 candentgroup.ca

📍 Unit 201H . 70 East Beaver Creek Rd
Richmond Hill . ON L4B 1L3 . Canada

