

# V-Blue

V e r s a t i l e   E n d o   F i l e

**NIC<sup>®</sup>**

For Better Endo

  
**CAN** ) **ENT**  
GROUP INC.

# V-Blue


V e r s a t i l e   E n d o   F i l e



Radiopaque stopper




Reciprocating

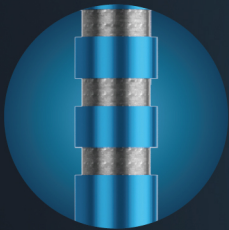

Precise depth marking



S Shape

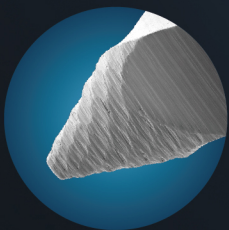

Safe guiding tip



# Features

- MaxTech 6.0 NiTi alloy
- S-shape cross section, strong centering ability
- Variable taper design
- 2 sharp cutting edge to create efficient cutting
- Reciprocating motion neutralize applied force to bring more safety
- Cover most morphology canals
- One file system



|              |               |               |               |
|--------------|---------------|---------------|---------------|
| V-Blue       | R25           | R40           | R50           |
| Size         | 25/.08        | 40/.06        | 50/.05        |
| Speed(RPM)   | Reciprocating | Reciprocating | Reciprocating |
| Torque(N.CM) | Reciprocating | Reciprocating | Reciprocating |
| GP           | R25           | R40           | R50           |

# Sequence

## Function 1: Shape and Finish Canals

- Establish the straightline access to canal orifice
- Choose the proper size V Blue according to X-ray film
- Shaping the canal by reciprocating 2/3. Better irrigate it every 3-4 pecks.
- Repeating shaping until it reaches working length

SSK #15



V-Blue



## Function 2: Canal Retreatment

- Remove obturation materials in canal coronal 1/3 with gates drills or any equivalent function device
- Rotate R25 by brushing the canal until working length

Gates drills





V-Blue



## References

| Size     | 21 mm       | 25 mm       | 31 mm       |
|----------|-------------|-------------|-------------|
| 25       | BFV 021 025 | BFV 025 025 | BFV 031 025 |
| 40       | BFV 021 040 | BFV 025 040 | BFV 031 040 |
| 50       | BFV 021 050 | BFV 025 050 | BFV 031 050 |
| 25/40/50 | BFV 021 A01 | BFV 025 A01 | BFV 031 A01 |

## Clinic Cases

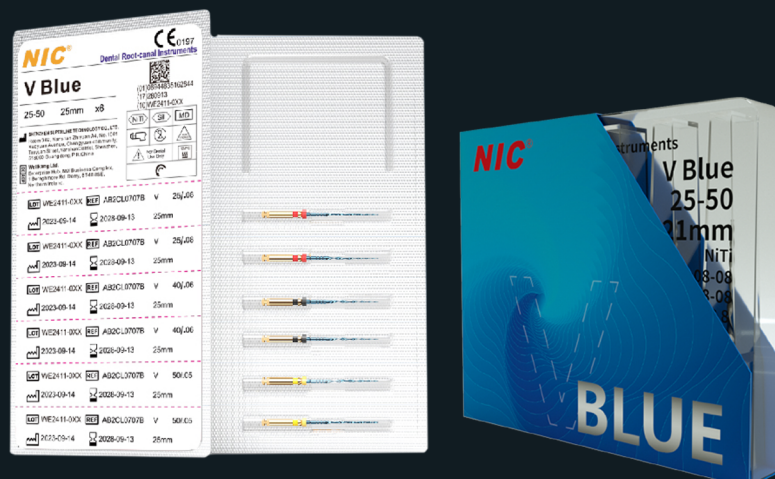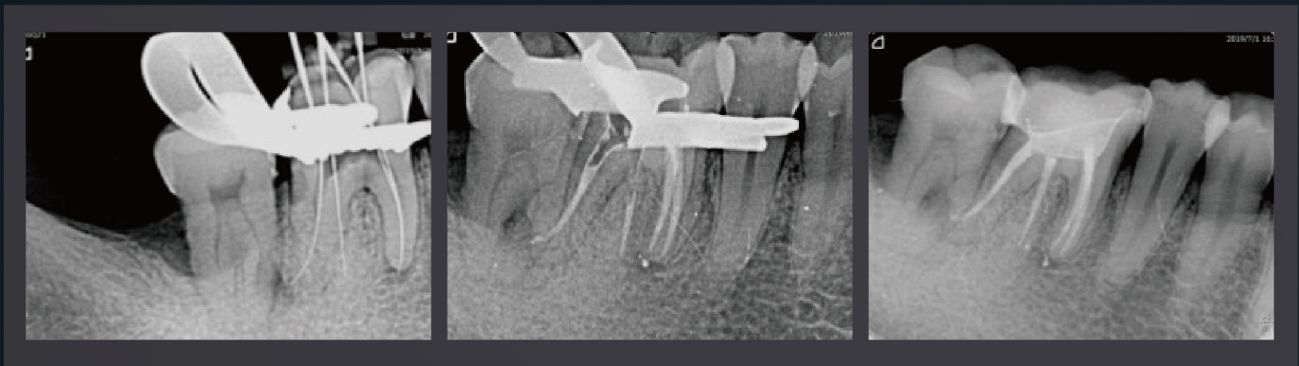

Available in sterile pack.



☎ 437-347-7791  
✉ [info@candentgroup.ca](mailto:info@candentgroup.ca)  
🌐 [candentgroup.ca](http://candentgroup.ca)

📍 Unit 201H . 70 East Beaver Creek Rd  
Richmond Hill . ON L4B 1L3 . Canada

