



CAN **DENT**
GROUP INC.

W + File

O n e F i l e E n d o

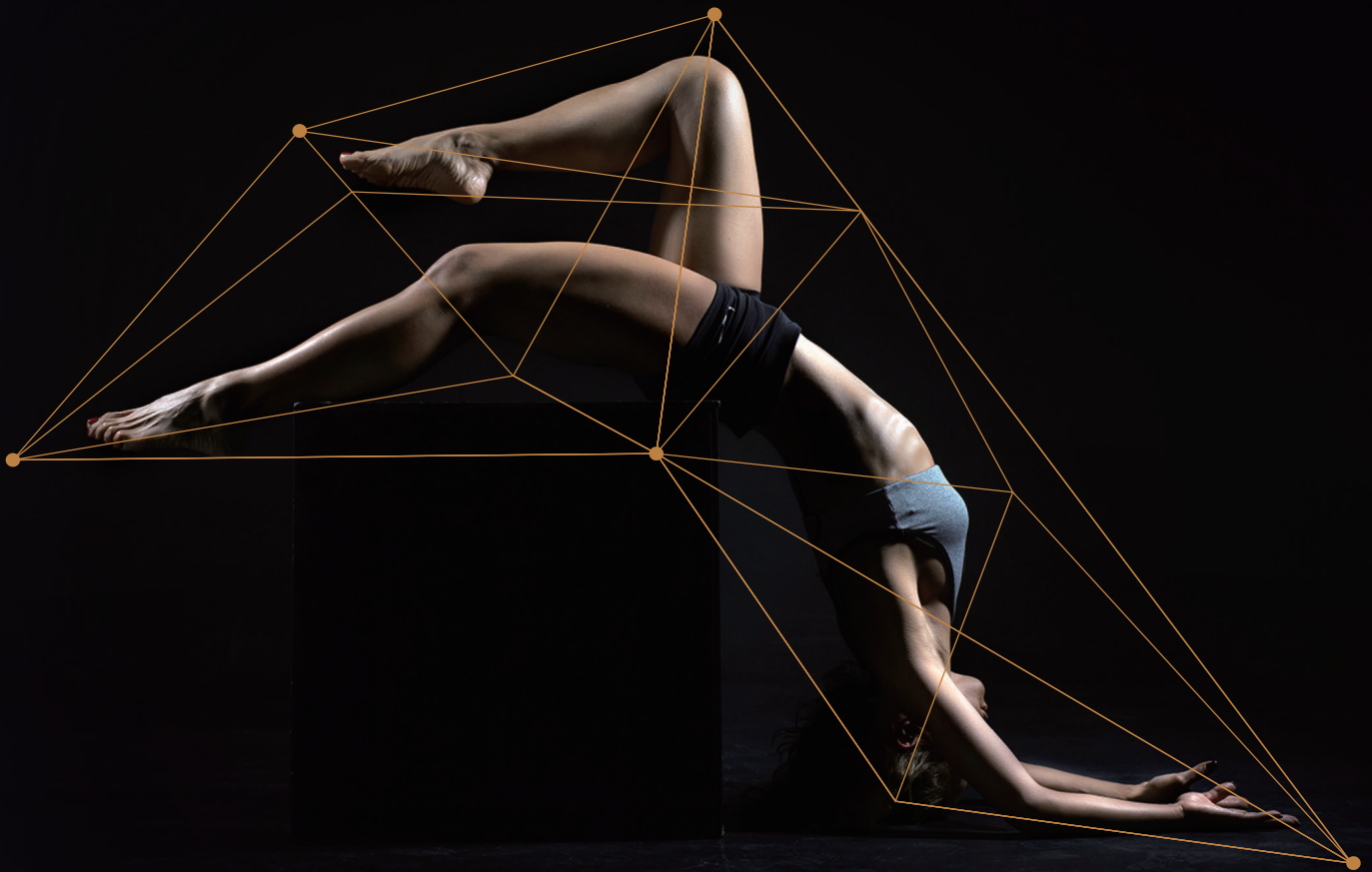

NIC[®]

For Better Endo

W+File

One File Endo




Radiopaque stopper




Precise depth marking




Reciprocating

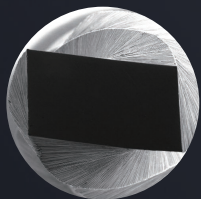

Parallelogram

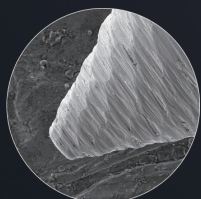

Safe guiding tip

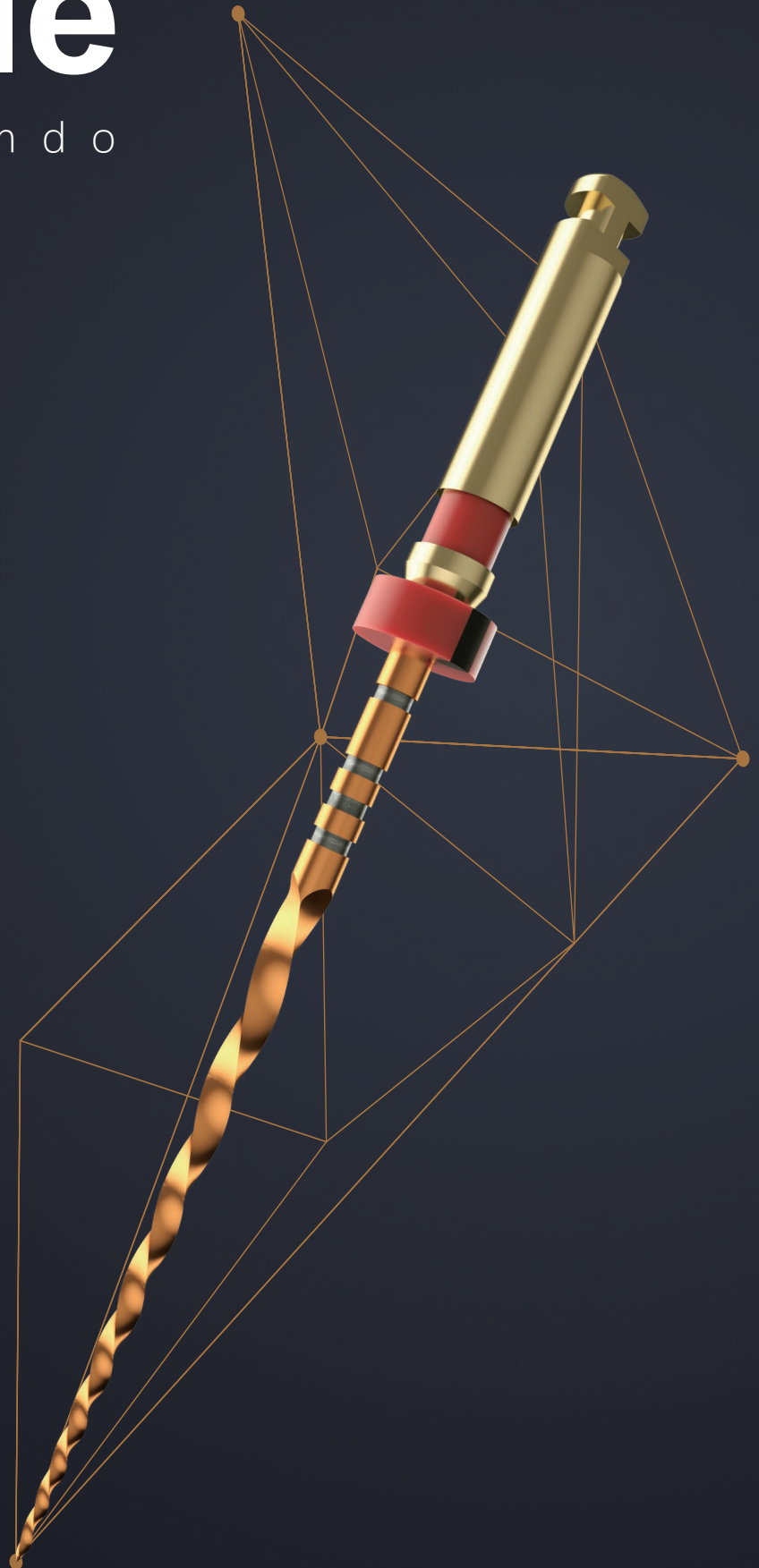

Features

- MaxTech 3.0 NiTi alloy
- Parallelogram cross section
- Variable regressive taper design to have accurate shaping
- Strong cutting ability shorten cutting time
- More space for debris removal
- Progressive pitch reduce screw-in
- One file system

				
W+File	Small	Primary	Medium	Large
Size	20/.07	25/.07	35/.06	45/.05
Speed(RPM)	Reciprocating	Reciprocating	Reciprocating	Reciprocating
Torque(N.CM)	Reciprocating	Reciprocating	Reciprocating	Reciprocating
GP	Small	Primary	Medium	Large

★ If Speed is adjustable under reciprocating motion,V file is recommended to set at:300-350RPM.

Sequence

- Scout the canal with hand use file
- Create glide path by any path function file
- Choose proper size shaping files
- Shaping the canal with the chose W+File until it reaches working length
- After shaping and finishing is done, then choose the matched GP points for obturation.

SSK #10



T+File



Primary



◆ Small 20/.07 for narrow canals

◆ Primary 25/.07 for most cases

◆ Medium 35/.06 & Large 45/.05
for wide canals



Small
20/.07



Medium
35/.06







Large
45/.05

References

Size	21 mm	25 mm	31 mm
20	AFW 021 020	AFW 025 020	AFW 031 020
25	AFW 021 025	AFW 025 025	AFW 031 025
35	AFW 021 035	AFW 025 035	AFW 031 035
45	AFW 021 045	AFW 025 045	AFW 031 045
20/25/35/45	AFW 021 A01	AFW 025 A01	AFW 031 A01

Clinic Cases

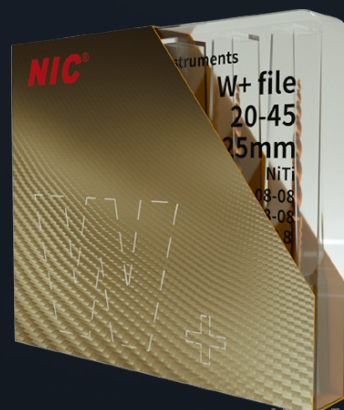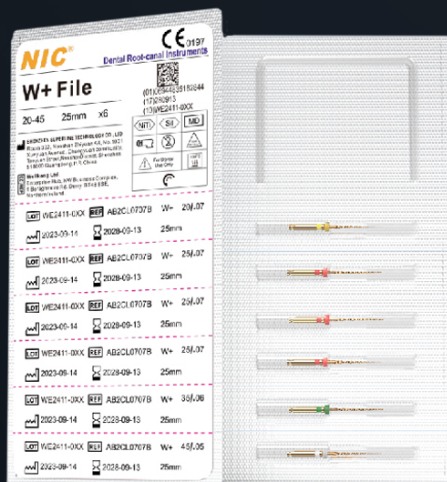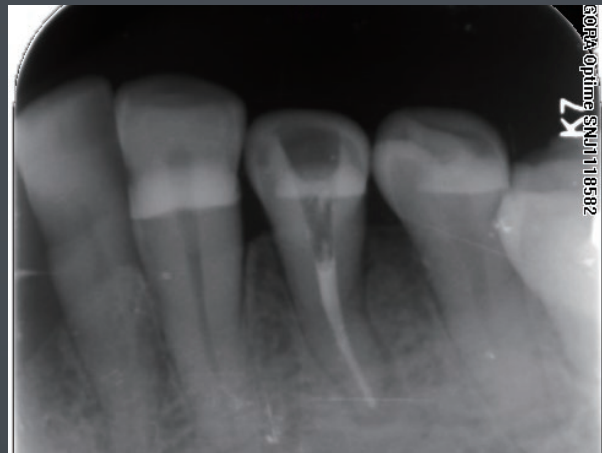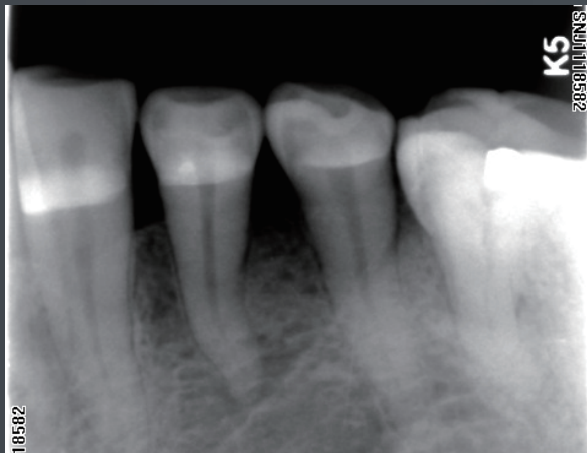

Available in Sterile pack.



☎ 437-347-7791

✉ info@candentgroup.ca

🌐 candentgroup.ca

📍 Unit 201H . 70 East Beaver Creek Rd
Richmond Hill . ON L4B 1L3 . Canada

